



clic  **pen**[®]
AMPOULE OPENER

SEPHA



CLIC-OPEN is a novel safety gadget for medics and patients who use ampoules and is ideal for pharmaceutical brand promotion.

FEATURES

- Easy to use
- Reduces injury
- Hygienic
- Sterilizable
- Reduces contamination

HOW DOES IT WORK?

Ampoules are usually opened by hand, using a file or with protective gauze. As the head of the ampoule is snapped off, the glass often breaks leaving dangerous, jagged edges which can cut the user's fingers. Slivers of glass or gauze fibers can also contaminate the liquid drug.

CLIC-OPEN reduces these hazards by using a tiny tungsten carbide cutter to pre-score the glass. The ampoule head can then be snapped open without jagged edges.

AN EFFECTIVE MARKETING TOOL

CLIC-OPEN is an instant hit with nurses, doctors, paramedics and patients. Every time CLIC-OPEN is used, the colour logo printed on the side encourages brand re-enforcement and customer loyalty.


CLIC-OPEN is not only a novelty item but also a useful safety gadget, making it the ideal ethical promotional gift for the medical industry. It can be given away either on its own or as part of an injectable's kit.

INSTRUCTIONS

1. Insert the ampoule head into the jaws of the CLIC-OPEN. Holding the ampoule steady, press only the upper jaw with the cutter onto the ampoule neck. Turn lightly to score.
2. Without squeezing the jaws together, but using the lower jaw only as an anchor, break open the ampoule at the score line.
3. Push in the plunger to eject the head of the ampoule into a sharps' waster container.



1 Score ampoule lightly with about 1/4 turn of CLIC-OPEN

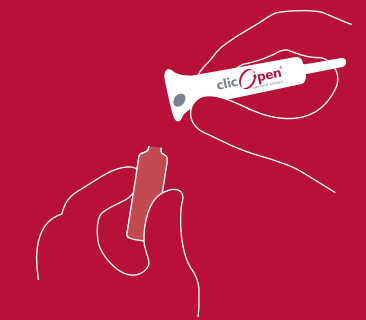

2 Snap with cutter opposite thumb



3 Eject head of ampoule



PRODUCT SPECIFICATION

- Suitable for ampoules up to size 20ml DIN
- Each CLIC-OPEN opens up to 200 ampoules
- **Dimensions**
Length: 64mm (2.5 inches)
Length with plunger extended: 102mm (4 inches)
Width: 25.4mm (1 inch) diameter at widest point
- Each CLIC-OPEN is supplied with a multi-lingual instruction leaflet (English, French, Spanish and German)

PRINTING

- Minimum order quantity for printing is 500 units
- Please contact Sepha for costs, stating number of colours
- Print area available: 12.7mm x 38mm (0.5 x 1.5 inches)
- Artwork specification: please supply film positive, emulsion side down

QUANTITY & PRICES

Product is available in the following quantities. Please contact Sepha for prices, or go to www.clicopen.com

QUANTITY	PACKAGING OPTIONS
5 - 49	Packed individually
50 - 499	Bundled in 100's
500 - 4,999	Bundled in 100's or Packed individually
5000 - 19,999	Bundled in 100's or Packed individually
20,000 +	Bundled in 100's

TERMS

Minimum order: Quantity 5 units
Delivery FOB (Belfast): 4-6 weeks from receipt of order, down payment and clearance of all details
Payment:
 Orders for 500 units or less: With Order
 Orders for over 500 units: 50% With Order
 50% on receipt of goods



INQUIRY FORM

For further information, please complete this section, detach along marked cutting line and return by:
Fax: 44 (0)28 9182 4255 or **Post:** SEPHA, Unit 25 Carrowreagh Business Park, Carrowreagh Road, Dundonald, Belfast BT16 1QQ, United Kingdom

Your Company Details:

NAME _____

JOB TITLE _____

COMPANY/ORGANISATION _____

ADDRESS _____

POSTCODE _____ COUNTRY _____

TEL FAX _____

EMAIL _____ WEBSITE _____

Someone else you would like us to send information to:

NAME _____

JOB TITLE _____

COMPANY/ORGANISATION _____

ADDRESS _____

POSTCODE _____ COUNTRY _____

TEL FAX _____

EMAIL _____ WEBSITE _____

What is your company's business activity?

HOSPITAL / CLINIC PHARMACEUTICAL COMPANY

ADVERTISING / PROMOTIONAL AGENCY

OTHER (PLEASE STATE): _____

Tick boxes relevant to your inquiry:

PLEASE SEND ME A FREE CLIC-OPEN SAMPLE

PLEASE CALL ME TO DISCUSS CLIC-OPEN FURTHER

I AM INTERESTED IN CLIC-OPEN FOR:

My own use

As a promotional item for my company

As a promotional item for a client

I would like a quotation for _____ units

SEPHA
working together

Unit 25 Carrowreagh Business Park, Carrowreagh Road,
Dundonald, Belfast BT16 1QQ, United Kingdom

T: +44 (0)28 9048 4848

E: sales@sepha.com

F: +44 (0)28 9048 0890

W: www.sepha.com

clic  **pen**[®]
AMPOULE OPENER